



60RBC200



YOUR
DIAGNOSTIC
COMPANION



Simple walk away operation In veterinary practices

- LCD touch screen
- Multi Languages
- Built in Plate Warmer
- Trusted Biogal Accuracy

- Protective Cover
- Magnetic Comb and paper strip holder
- Digital Instructions For Use
- Factory Calibrated






Working with RoboComb+ is simple

Select the desired test





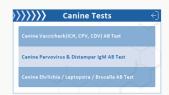

- Bring the plate to room temperature by: 2
 - 1. Leaving out of refrigeration for 60 minutes or
 - 2. Using the built-in plate warmer for 10 minutes
- Load the sample into the kit's plate and place the plate on the plate base of the RoboComb+
- Load the tooth, teeth or full comb on the comb holder

1-12 teeth can be tested simultaneously



- Let the test run
- Read the results







Designed for Biogal's antibody tests

COMPATIBLE WITH BIOGAL'S TESTS:



Validate your Patient's protection

Scan OR





- Infectious Hepatitis
- Parvovirus
- Distemper



- Panleukopenia
- Herpesvirus
- Calicivirus

ImmunoComb[®]

Antibody levels assessment

Scan QR





- Parvovirus & Distemper IgM
- Ehrlichia
- Leptospira
- Brucella



- Toxoplasma & Chlamydophila
- Coronavirus (FCoV)



- Brucella
- Leptospira



- IBD-ND-IB
- Mycoplasma Mg-Ms-Mm
- Neospora





Automate your test reading

CombCam

Objectively reads VacciCheck / ImmunoComb Antibody test results







